

IMPLANT SYSTEM

Neobiotech

IMPLANT SYSTEM

**Neo
Biotech**
Satisfaction to Dentists

“Türkiye Distribütörü”


AC DENTAL®
www.acdental.com.tr

The image features three wireframe trays of surgical instruments, each containing various tools like scalpels, forceps, and probes. The trays are arranged in a triangular pattern against a dark, smoky background. The top tray is at the top right, the middle one is in the center left, and the bottom one is at the bottom right. The text is overlaid on the middle tray.

NeoBiotech
SURGICAL KIT
PROSTHETIC KIT

Neo Master Kit


This kit can be used to perform implant surgery for IS-III active, IS-II active, IT-III active

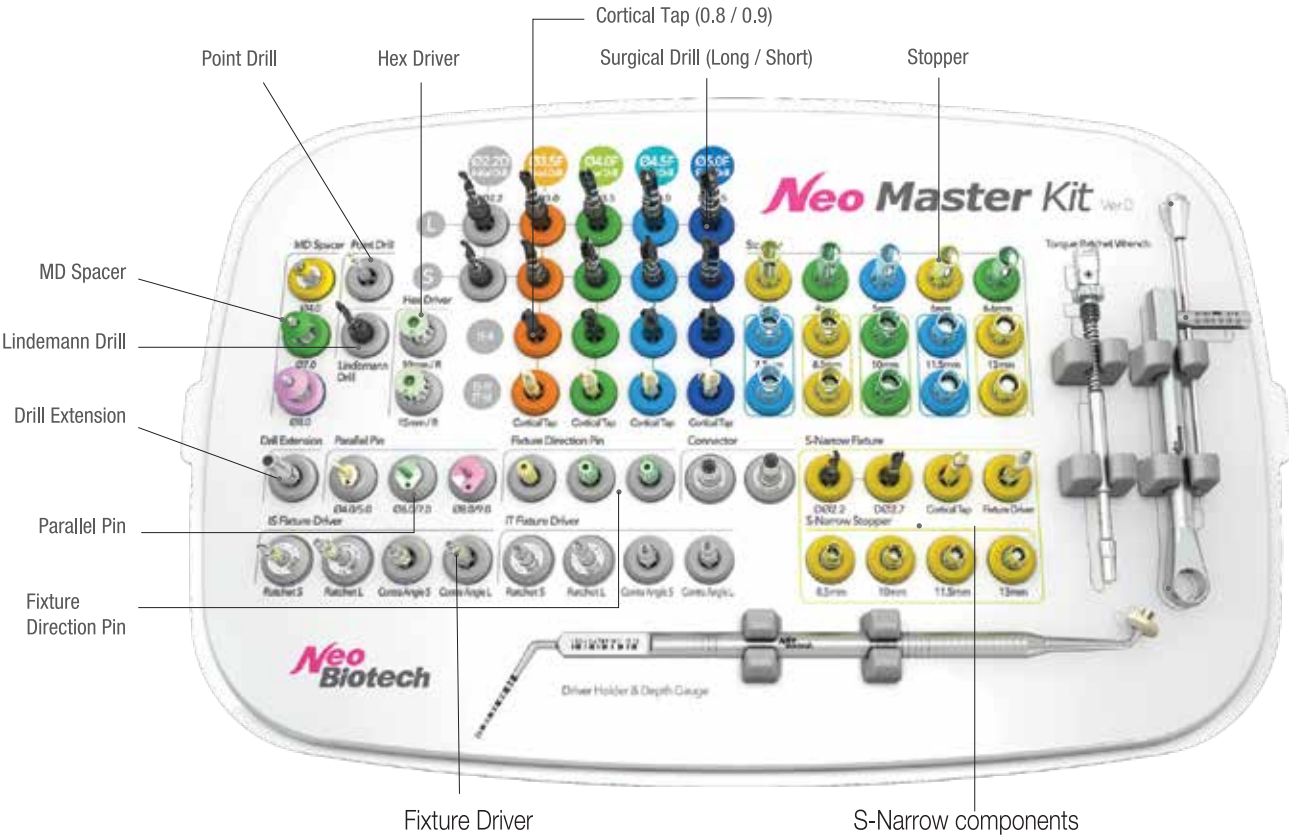

kit	NEOMKIT
Kit Tray Top	NEOMKITTRAY1
Kit Tray Middle	NEOMKITTRAY2
Kit Tray Bottom	NEOMKITTRAY3

Neo Core Kit

Simple Surgical Kit for IS-II active / IS-III active

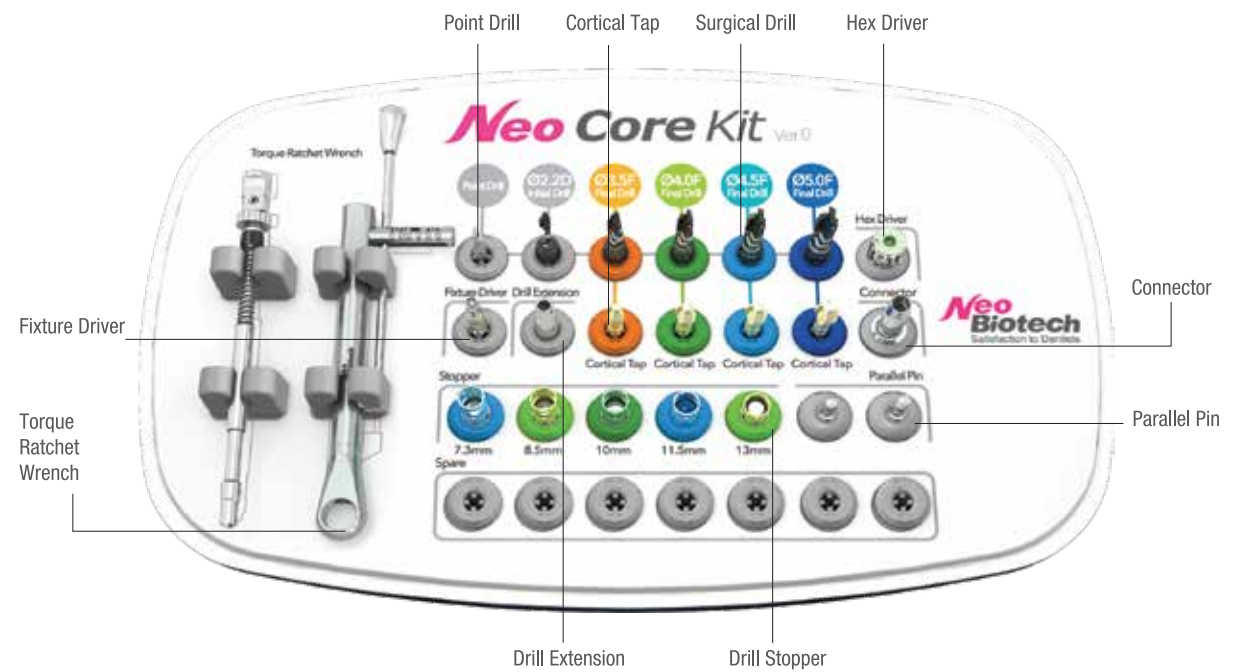


Kit	NEOCKIT
-----	---------

S-Narrow Kit

This kit can be used to perform implant surgery for IS-III active(Ø3.2).


• All products are individual purchases.

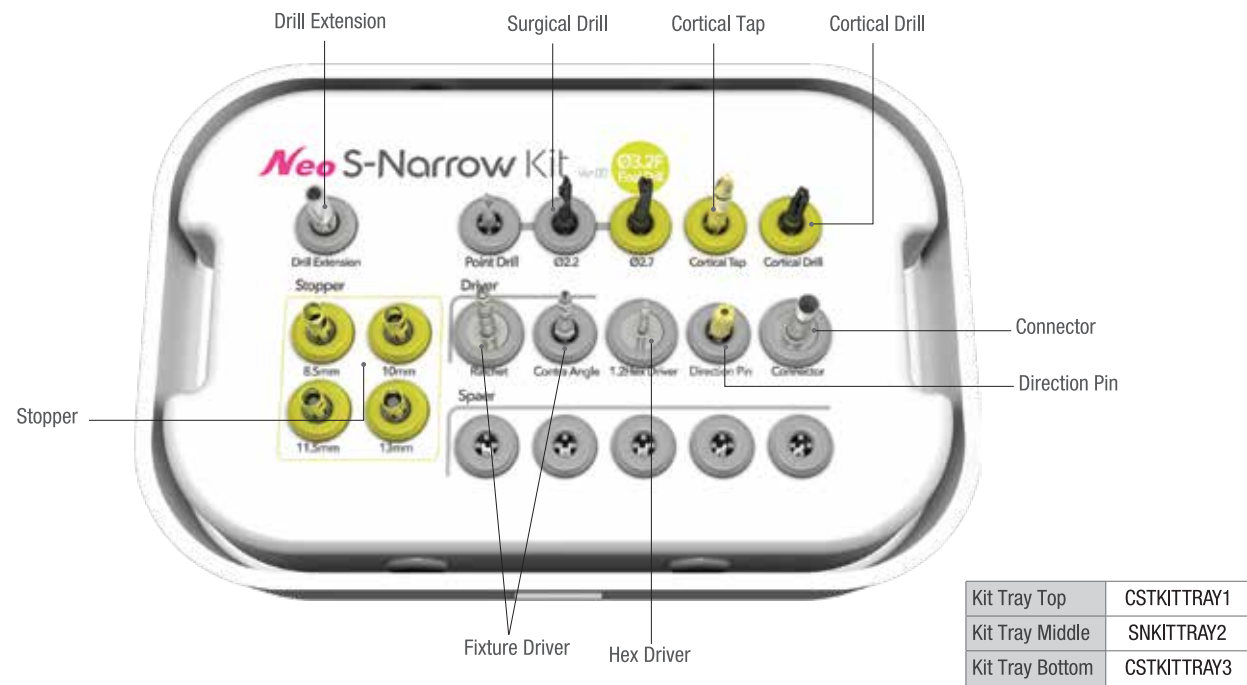

Kit Tray Top	CSTKITTRAY1
Kit Tray Middle	SNKITTRAY2
Kit Tray Bottom	CSTKITTRAY3

Neo NaviGuide Kit


This kit can be used to perform surgical guide implant surgery for IS-II active , IS-III active (Ø3.5~Ø5.0) For R/W sleeve , R/W C sleeve

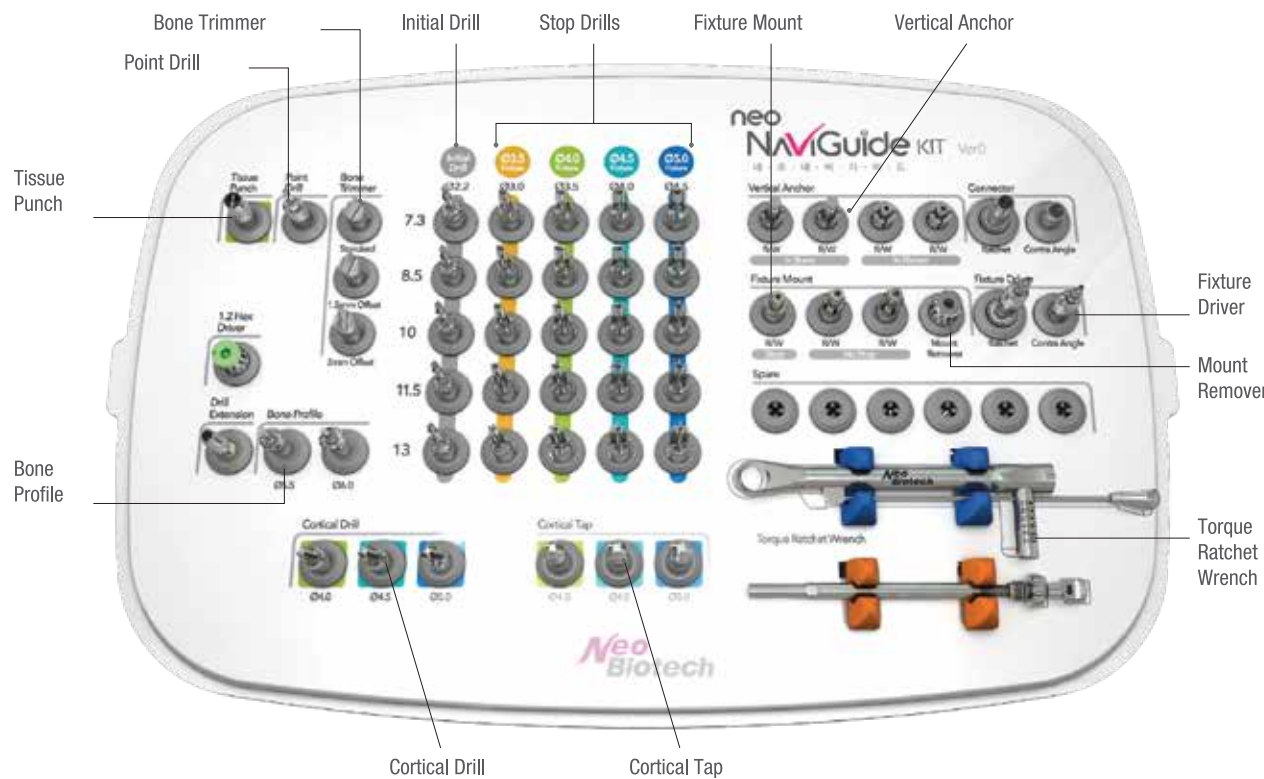


Product name	NEONGKIT
--------------	----------

Surgical Guide Sleeve


Narrow

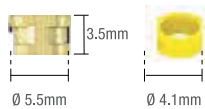

Regular/Wide

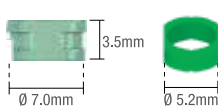

C-Sleeve Regular/Wide

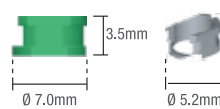


Diameter	5,5	7	7
Anti-rotational distance	4,1	5,2	5,2
Height	3,5	3,5	3,5
Product name	SGS55N	SGS70RW	CSGS70RW




Ridge Wider Kit

The product is a dental procedural device developed for convenient and easy split and expansion with RidgeWider Kit including expander, safe disk, drill and chisel in narrow ridge case where implant placement is difficult.

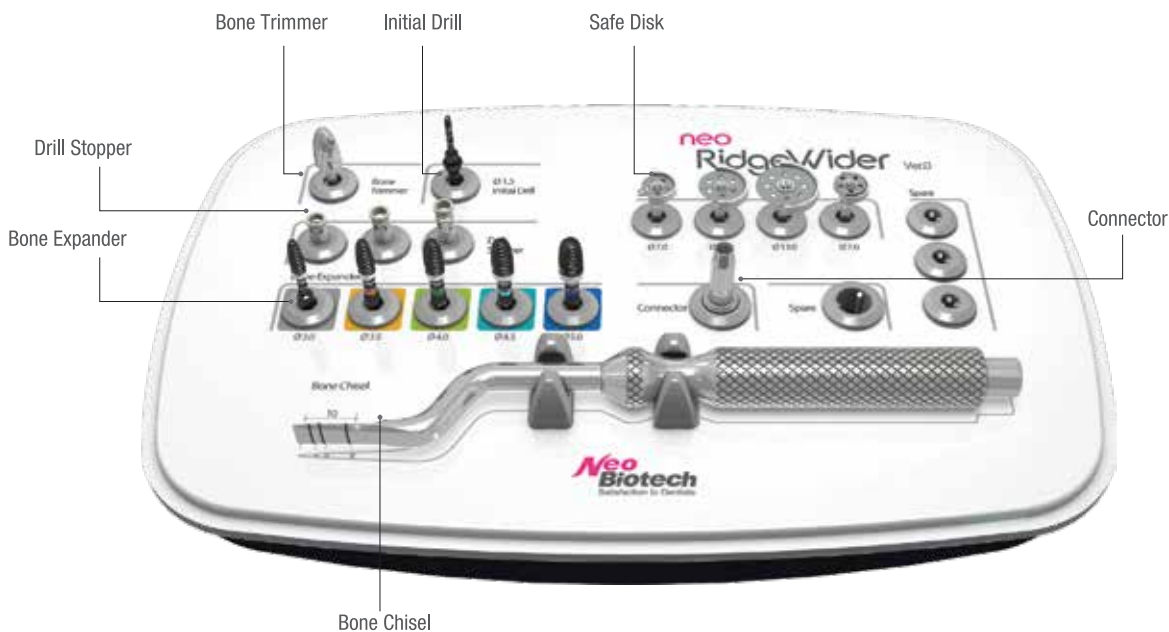


Kit	RWKIT
-----	-------

Neo Surgical Accessory Kit

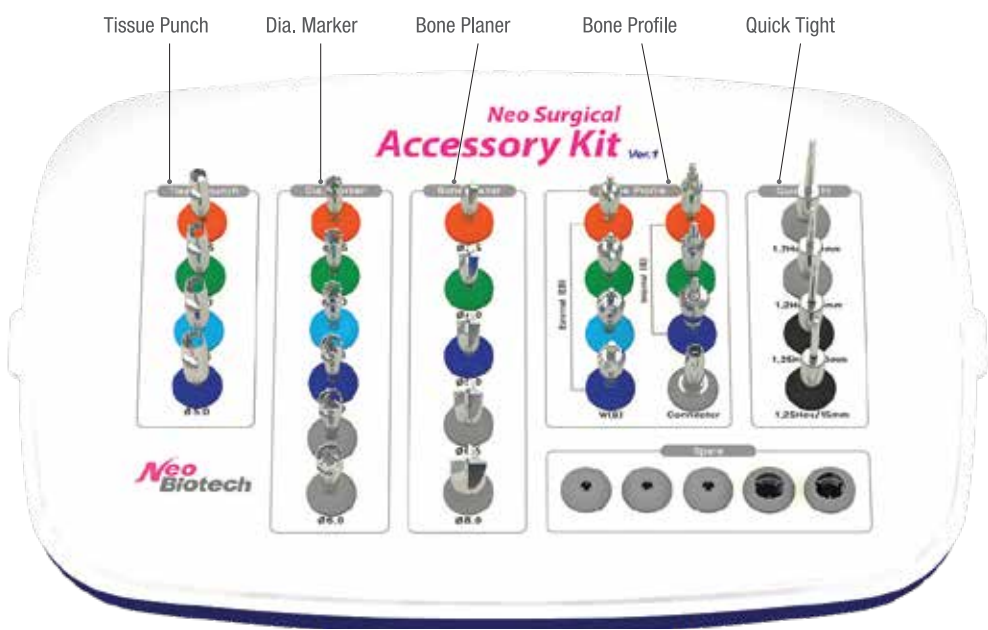


Kit	ACCKIT
-----	--------

Neo Prosthetic Kit

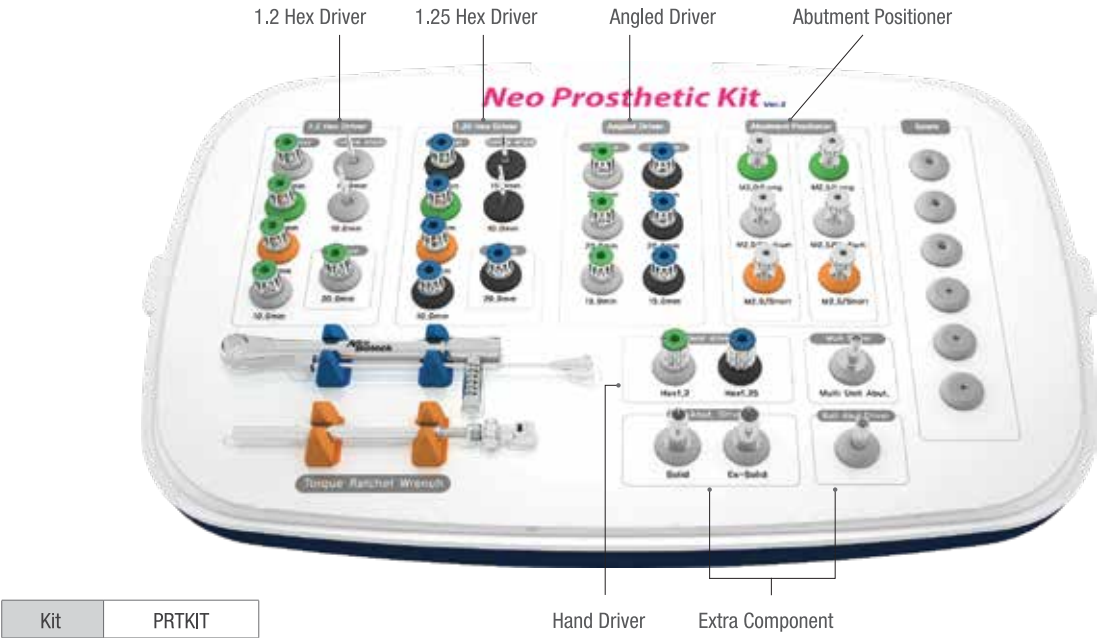


Neo Prosthetic Planning Kit

Neo Prosthetic Planning kit which helps users select a prosthesis by diameter, cuff and length when engaging with a fixture before choosing the size of prosthesis

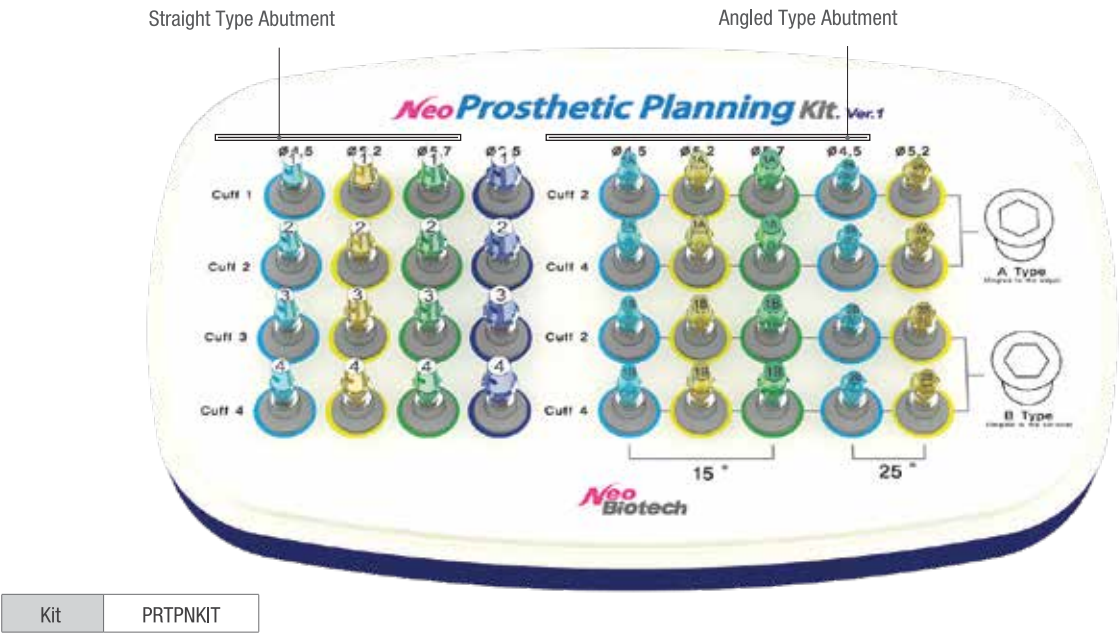

•Hex Driver

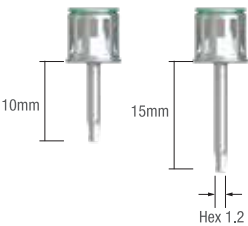

Length (mm)	Hex	Product name
10	1.2	HD1210S
15		HD1215S

• Torque Ratchet




Product name	TW60
--------------	------

A wireframe 3D rendering of a Neobiotech Solution System. The system consists of three main components: a large rectangular base unit with a grid of circular ports, a smaller rectangular unit with a grid of ports, and a central rectangular unit with a grid of ports. The components are interconnected by a network of tubes and connectors. The entire system is rendered in a wireframe style, showing the internal structure and connections. The background is dark with a subtle grid pattern.

Neobiotech
SOLUTION SYSTEM


SinusAll Kit

The SinusAll Kit combines NeoBiotech’s SCA (Sinus Crestal Approach) Kit, SLA (Sinus Lateral Approach) Kit, and the Kit into a single kit. The product allows maxillary sinus implant procedures to take place using a single kit, rather than to prepare a sinus kit and a surgical kit separately.

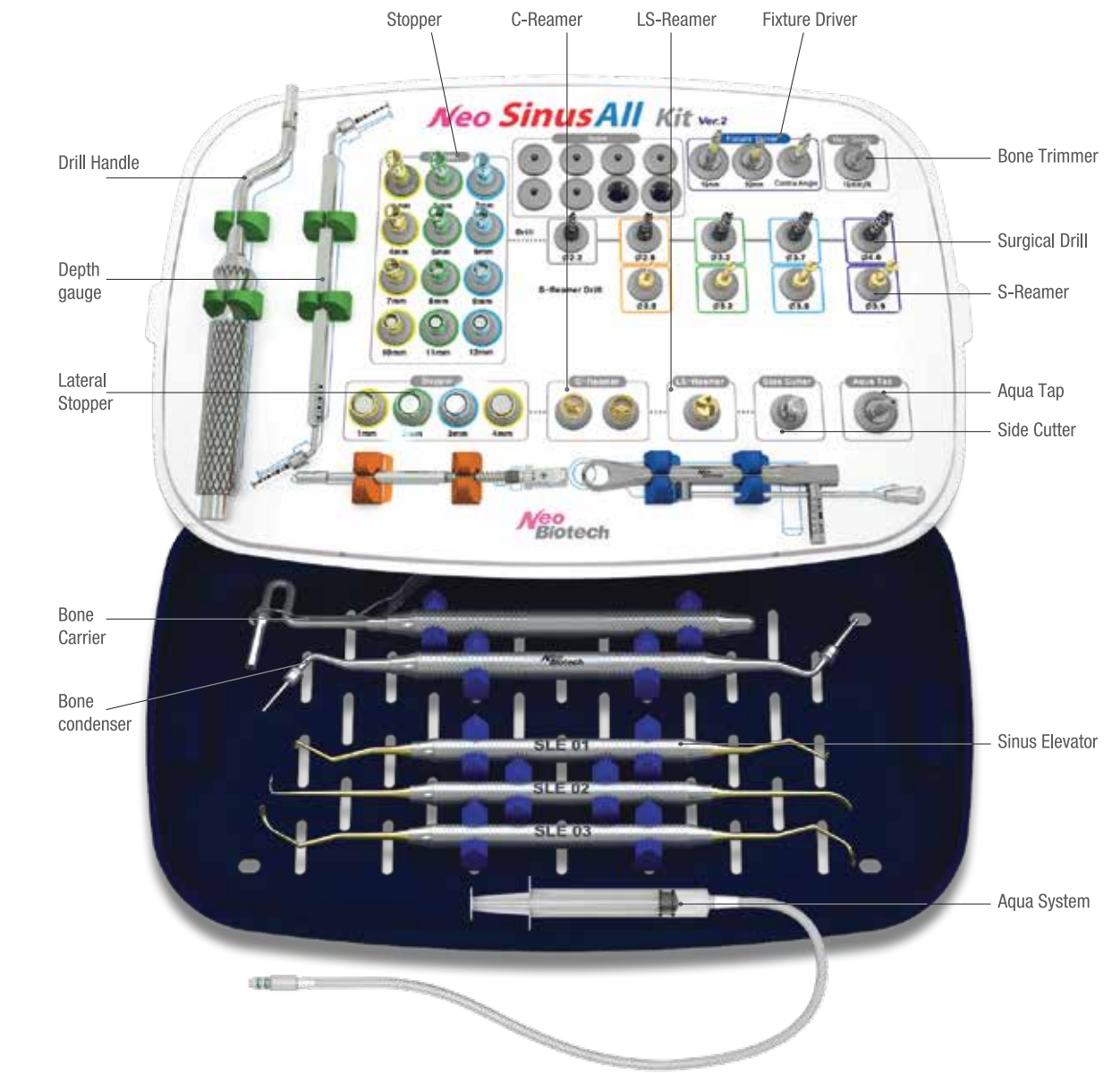


SCA Kit (Sinus Crestal Approach Kit)

SCA Kit is a drilling tool on the inferior cortical bone without tearing of membrane and make a perforation on the inferior cortical bone with out malleting osteotome technique Therefore User can get CMI Fixation. (Initial Fixation in Crest cortical bone, Middle cancellous bone and Inferior cortical bone)

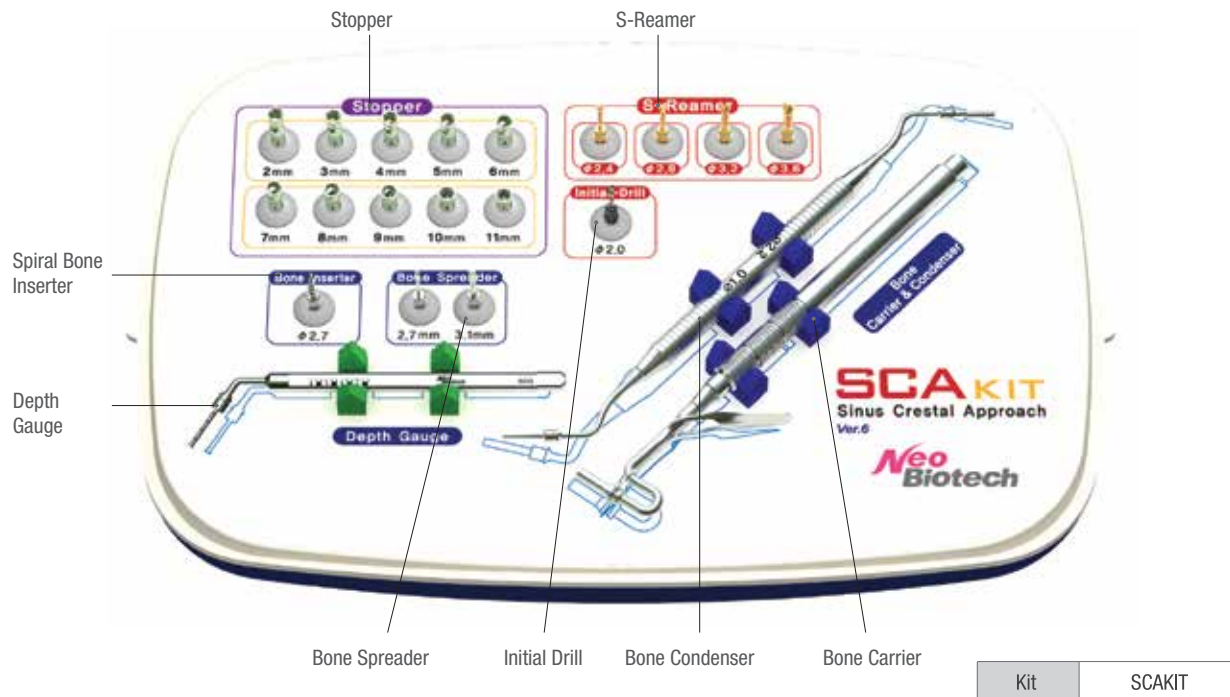


SLA Kit (Sinus Lateral Approach Kit)

Lateral window technique is the formation of the access to the maxillary sinus through its lateral wall. The process tool is simpler and much comfortable than any other techniques. Above all, this surgical kit provides ultimate solutions for the cases of only having with a thin residual bone height which is difficult to approach with a crestal technique, membrane rupture by a sinus crestal approach or placing a multiple implants. In addition, SLA reamers can be easily adapted for the case of formation of minimum flap size on the inferior wall.

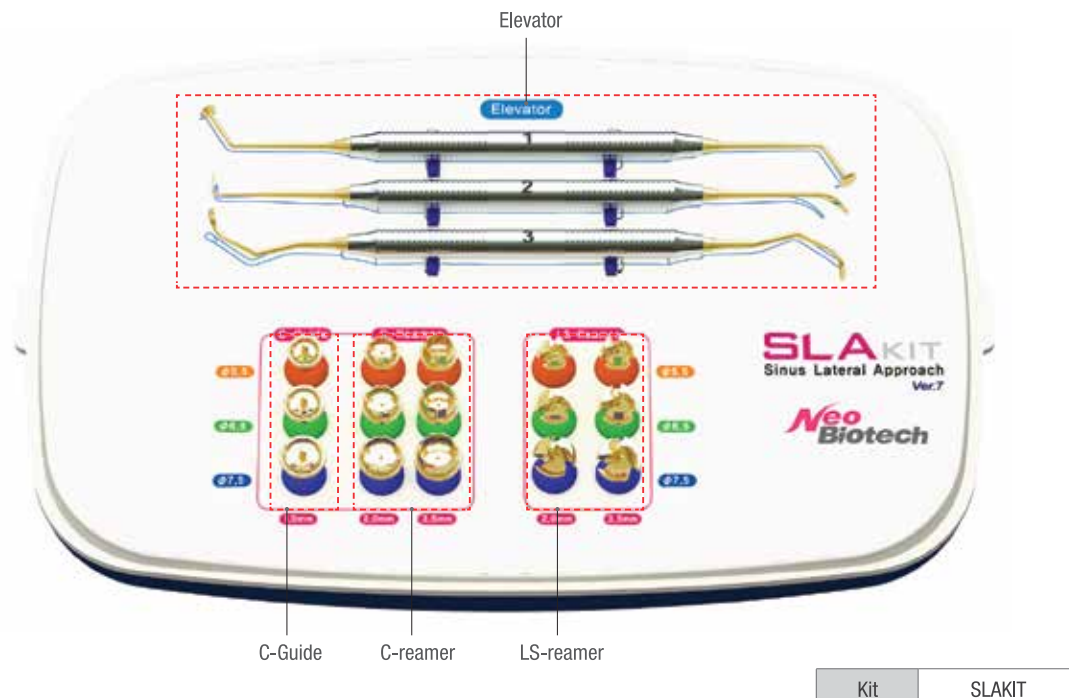


FR Kit (Fixture Remover Kit)

This product is a surgical tool designed to remove implants that were stopped during the implantation due to torque or implants whose surrounding bones have been damaged. After removing the implant, a new implant with the same diameter can be immediately implanted.

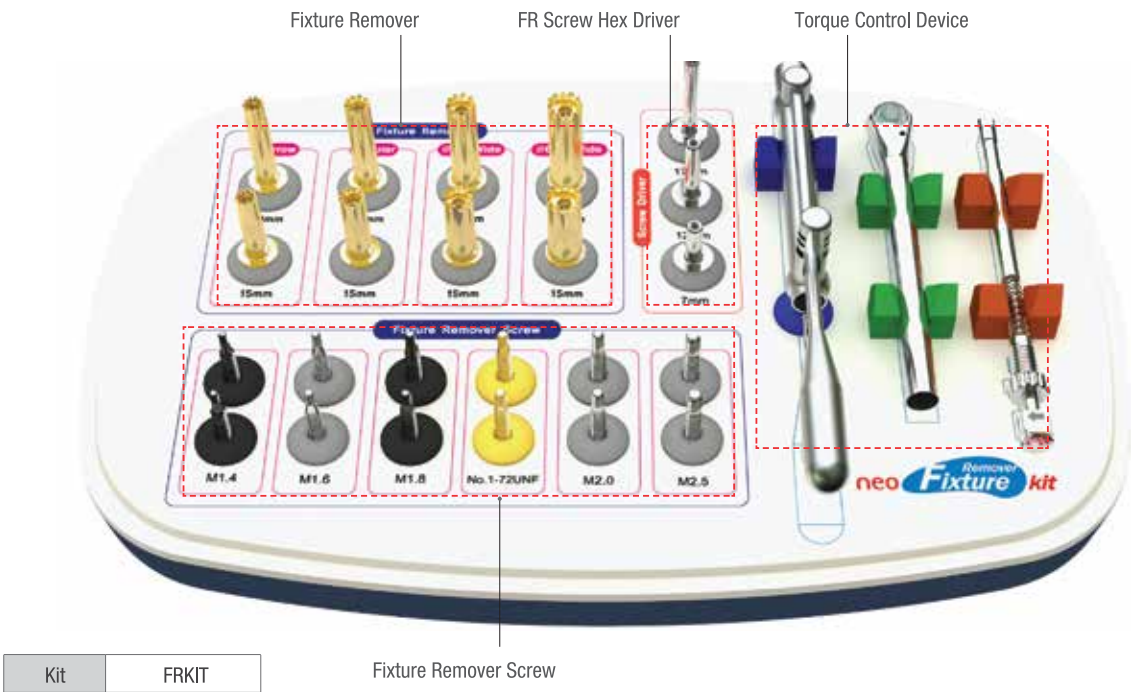

SR Kit (Screw Remover Kit)


This product is a surgical tool developed in order safely and quickly to remove screws that have become fractured inside an implanted Fixture for various reasons. After removing the Screw, a new Abutment may be connected to the Fixture.

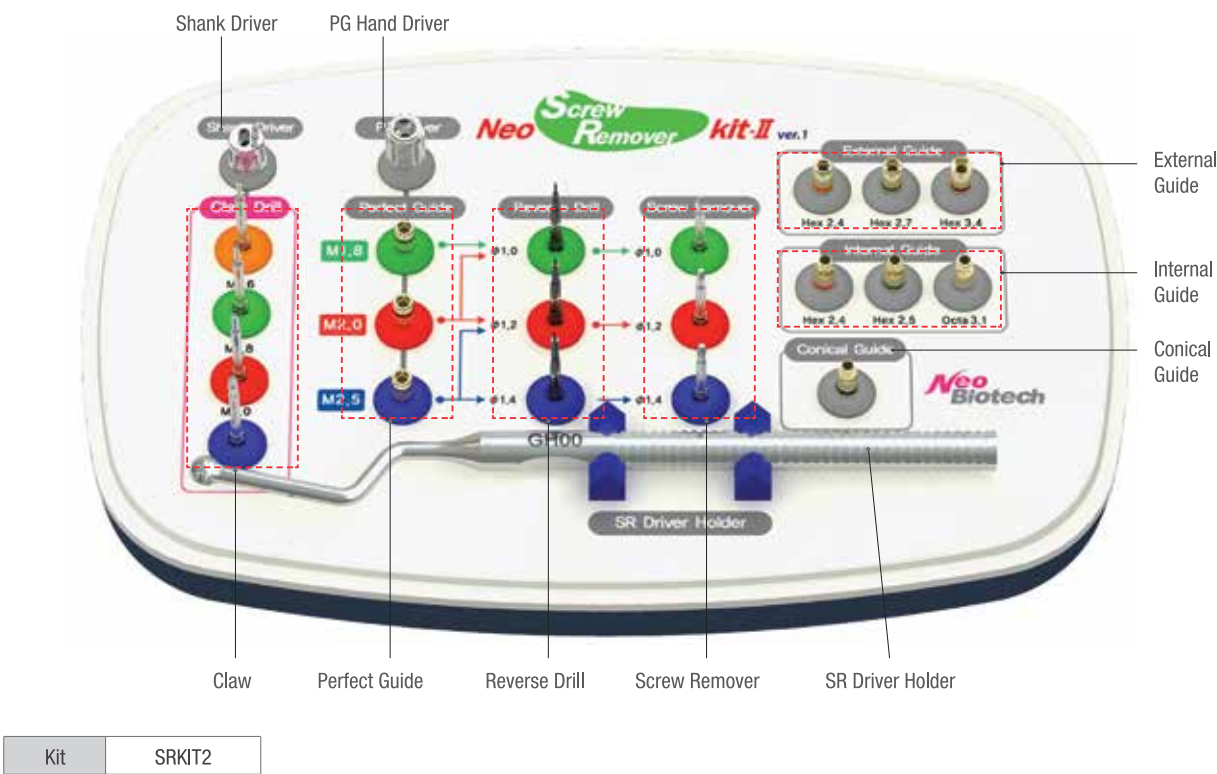

ACM (Auto Chip Maker)

ACM is an instrument to collect autograft bone from mandibular molars or incisor for GBR in defected area.

• Auto Chip Maker

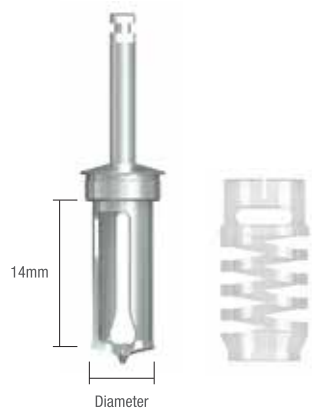

Diameter (Ø)	Product name
Ø 4.0	NACM40SETS
Ø 4.5	NACM45SETS
Ø 5.0	NACM50SETS
Ø 6.0	NACM60SETS

※ ACM Stopper 7ea Included

• Auto Chip Maker Technique Guide



Stopper

- Connect the stopper to the sterilized ACM drill. .



Drilling

- To prevent slipping of the ACM drill, slightly tilt the drill left and right until the drill stably positions itself. (80rpm / Non-irrigation recommended)
- In order to collect high quality autograft bone, drill only 3~4mm deep from on site, and move on to a different site to continue drilling.



Locate the ACM

- Select the most suitable place to collect autograft bone.
- Position the ACM Drill where the autograft bone will be collected.




EZ-GBR Kit

EZ-GBR Kit is a surgical tool used for oral bone grafting onto areas in which there are defects in autogenous bones. It consists of surgical tool components that can affix membranes.




CTi-mem


New Concept of Titanium Membrane, Customized Ti Membrane

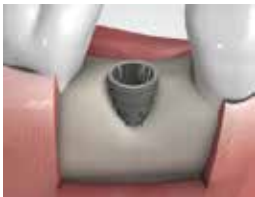


• Bone defect in Buccal area

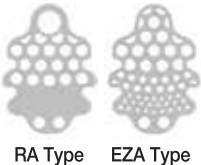


• Bone defect in Buccal, Proximal area

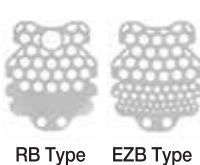

• Bone defect in Buccal Proximal, Lingual area

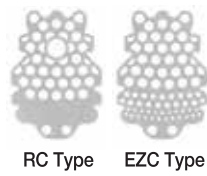


• Bone defect large area

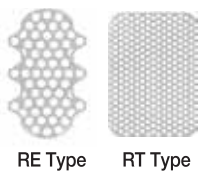

• Tent Screw Type



R Type CTi-mem

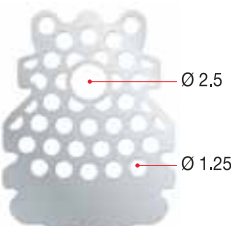


Thickness : 0.12mm

• EZ-Fixing Screw Type



EZ Type CTi-mem

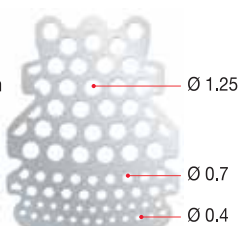

Thickness : 0.12mm

AnyCheck

AnyCheck is an implant stability meter that measures the stiffness of the alveolar bone-implant interface through a tapping-motion.

AnyCheck reduces the risk of implant failure by indicating the most appropriate time for connecting prosthetics.



Measured value is similar to ISQ scale: above 60 – bridge recommended; above 70 – single case loading possible.

Displayed in red – 30~59 ; orange – 60~64 ; green – 65~85.

Measured values (IST scale) differ depending on the height of the healing abutment. As the standard height of the healing abutment is 4mm, calculate the correct IST value according to the following chart:

Standard +3	7mm height Healing Abutment	Displayed IST value +6
Standard +2	6mm height Healing Abutment	Displayed IST value +4
Standard +1	5mm height Healing Abutment	Displayed IST value +2
Standard	4mm height Healing Abutment	Displayed IST value
Standard - 1	3mm height Healing Abutment	Displayed IST value - 2
Standard - 2	2mm height Healing Abutment	Displayed IST value - 4
Standard - 3	1mm height Healing Abutment	Displayed IST value - 6

Peri-implantitis Solution

• R-Brush



Product Code	
Regular (Ø5.6)	ROUND BRUSH (R)
Wide (Ø6.6)	ROUND BRUSH (W)

• i-Brush 1



Product name	IBRUSH1
--------------	---------

• i-Brush 2



Product name	IBRUSH 2
--------------	----------



Neobiotech Implant System